

The Steraffirm Steam Class 6 Emulating Indicators are cycle verification indicators that confirm the presence or absence of specific time and temperature parameters during a cycle, and integrate all the critical parameters of the steam sterilization cycle (temperature, saturated steam and exposure time).



Steraffirm® Steam Class 6 Emulating Indicators

FEATURES

Abrupt and distinctive total color change from yellow to blue/purple

Endpoint color standard is present on each indicator strip

No lead or other heavy metals

Indicator can be kept as a permanent record

Lot number on each indicator strip and on each box

Conforms to ISO 11140-1 (Class 6)

BENEFITS

Easy to interpret test results

Indicates when time and temperature parameters are achieved

No concern over exposure to lead or other heavy metals


Valuable documentation of the efficacy of the sterilization cycle

Enhances traceability

Meets critical standards as an indicator for monitoring specific cycle parameters


Emulating Indicators can be placed inside each pack or pouch during sterilization of supplies or equipment in a validated prevacuum cycle. The indicators must be checked by the end user (before using the processed item) to confirm that the specific sterilizing parameters have been attained within the pack. Class 6 emulating indicators can also be placed in test packs and processed in a cycle. The test packs should be opened immediately after the cycle to assess the success of the process. These emulating indicators can be used for routine monitoring of research or pharmaceutical steam sterilization cycles. The Steraffirm Class 6 Emulating Indicators are for life sciences applications only. They are not approved for use in healthcare facilities in the United States.

COLORS ARE AN APPROXIMATION




Steraffirm™ **STEAM** Emulating Indicator

RE ORDER No. PCC012 Class 6
For [STEAM] sterilizers working at 121°C / 250°F for 15 minutes
ACCEPT ONLY IF YELLOW CIRCLE TURNS BLUE/PURPLE




UNUSED




Steraffirm™ **STEAM** Emulating Indicator

RE ORDER No. PCC012 Class 6
For [STEAM] sterilizers working at 121°C / 250°F for 15 minutes
ACCEPT ONLY IF YELLOW CIRCLE TURNS BLUE/PURPLE



TYPICAL PASS

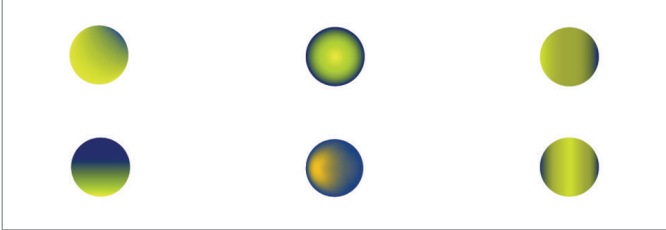

TYPICAL FAIL

PRODUCT	DESCRIPTION	QUANTITY
PCC012	Steraffirm Steam Class 6 Emulating Indicator 121°C/250°F for 15 minutes	100/box



STERIS CORPORATION
5960 Heisley Road
Mentor, OH 44060, USA
T +1 800 444 9009

www.STERISLifeSciences.com

STERIS Offices Worldwide

Asia Pacific	+65 (0) 68 41 7677
Canada	+1 800 444 9009
Finland	+358 9 25851
France	+33 (0) 2 38 84 85 40
Germany	+49 (0) 221 466120
Italy	+39 (0) 2 21303 424
Latin America	+1 800 444 9009
Spain	+34 (0) 916 585 920
United Kingdom	+44 (0) 1256 840400

STERIS Life Sciences | Science & Solutions for Life

For over 100 years, STERIS has been known as a global leader, trusted partner and solutions provider in the field of sterilization and contamination control. **Today**, STERIS continues building on this heritage by advancing the **Science** of sterilization, cleaning and infection control while offering **Solutions** that meet our Customers' needs and high standards. STERIS is dedicated to helping you enhance the **Life** of your patients and the life of your equipment. From Formulated Chemistries, to Capital Equipment, to Parts and Services, STERIS Life Sciences is **Science & Solutions for Life**.

© 2015 STERIS Corporation
All rights reserved

No part of this publication may be reproduced in any material form without the permission of STERIS Corp.